

INFORMATION BULLETIN

March 2026

NEW PERFORMA CG AIR THERAPEUTIC AIR CUSHIONS

INTRODUCTION

In September 2025, we proudly launched the **Performa CG Air Therapeutic Cushion**, the latest evolution of our CG Air technology. This new cushion maintains the same trusted performance our users and clinicians rely on, while introducing key improvements: a redesigned, softer cell structure, more compartment configuration options, and new breathable cover choices. These updates deliver even greater comfort, stability, and skin protection, helping us provide better outcomes for users and a stronger solution for our clinical partners.

OLD DESIGN VS NEW DESIGN

Previous CG Air Design



NEW Performa CG Air



NEW PERFORMA CG AIR CUSHIONS (cont'd)

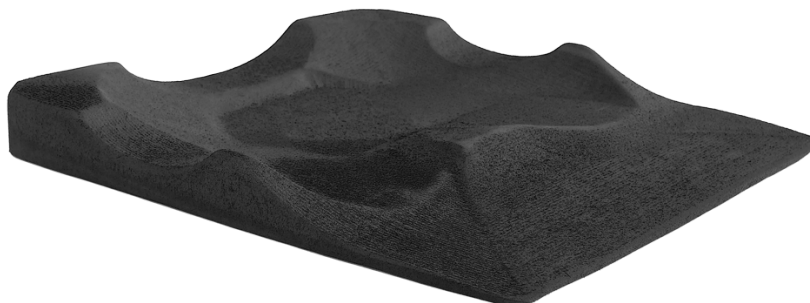

Component	Previous CG Air Design	NEW Performa CG Air
Cell Design and Material	<ul style="list-style-type: none"> • Medical-grade PVC material (non-allergic and anti-fungal) • Non-porous material resistant to discoloration and odors, with improved air retention. • High frequency welding and injection molding • Made from fire/flame retardant (CAL117-2013) 	<ul style="list-style-type: none"> • DINP-free and softer to the touch medical-grade PVC (non-allergic and anti-fungal) • Non-porous material resistant to discoloration and odors, with improved air retention • 15% thinner upper cell wall at the top for better immersion with consistent base support • High frequency welding and injection molding • Made from fire/flame retardant (CAL117-2013) • Resists sticking from dried nutritional supplement spills
Compartment Options and Cushion Configurations	<ul style="list-style-type: none"> • Dual valves allow for side-to-side positioning options • Custom dimensions and cell configurations, asymmetrical seat depth 	<ul style="list-style-type: none"> • Dual-compartment design: Left/right or front/rear configurations for personalized positioning and pressure management • Custom dimensions and cell configurations, asymmetrical seat depth
Cell Base	Maintains the same reliability and durability at the cell base	Maintains the same reliability and durability at the cell base
Covers	<ul style="list-style-type: none"> • Standard 2-way stretch • 3D mesh • Incontinence 	<ul style="list-style-type: none"> • New covers: improved breathability for better temperature and moisture control • Standard single-layer spacer knits (4-way stretch) • Optional 3D double-layer-spacer knits (4-way stretch) • Incontinence • NEW anti-slip technology (Patent Pending) option: GripZone – Prevents shear forces in the ischial and coccyx regions by placing a thin anti-slip polymer under the femurs and trochanters. <ul style="list-style-type: none"> ○ No reduction in the fabric's stretch or breathability ○ Stitch-free design
Clinical Benefits	<ul style="list-style-type: none"> • Superior pressure management: Stimulates circulation & reduces pressure points • Enhanced postural stability: Ischial well, reduced pelvic migration 	<ul style="list-style-type: none"> • Prevents posterior pelvic tilt • Unsurpassed stability, effective for foot propellers • Superior pressure management: Stimulates circulation & reduces pressure points • Enhanced postural stability: Ischial well, reduced pelvic migration

NEW PERFORMA CG AIR CUSHIONS (cont'd)

New Optional Equipment

Anatomic Contour Seat Insert

Increases pelvic support, immersion in sensitive areas, and lateral stability.




Convenient Pouch

Perfect for a cellphone and essentials, with easy to manage magnetic flap.



Transfer Trellis

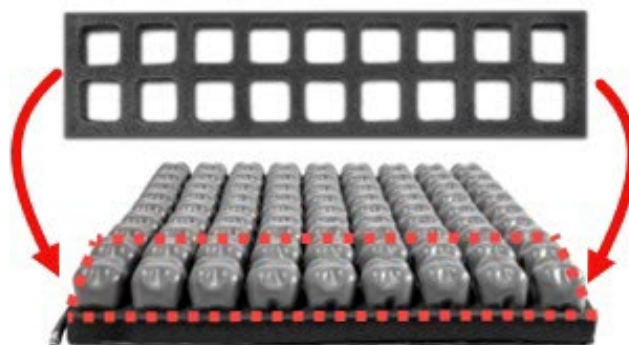Supports the user's weight during transfers in or out of the chair, protecting the vulnerable front cells. Creates enhanced front-end stability for foot propulsion.



NEW PERFORMA CG AIR CUSHIONS (cont'd)**OTHER CONSIDERATIONS****Availability, Pricing and Order Forms**

The Performa CG Air is now available, and pricing is listed in the updated Order Forms online. The Alltrack order forms have been updated with the new Performa CG Air information as well.

- US Customers: US Order Forms contain the revised HCPCS codes for the standard configurations.
- Ontario Customers: Ontario Order Forms contain the new ADP codes for all configurations.

Updates**Documents:**

All pertinent documents have been updated and are online.

- [Website](#)
- [Owner's Manual](#)
- [Order Forms USA](#)
- [Order Forms CA](#)
- [Order Forms ON](#)
- [Parts Manual](#)
- [Performa CG Air Flyer](#)

Videos:

Repair videos were created and are online under *Support / For Pros / [Instruction Guides - Videos](#)*.

- Repairing a tear at the base of the cell
- Repairing a detached valve
- Repairing a puncture from a sharp object