

It stimulates and helps the step pattern. It trains the user to maintain the balance of the trunk and pelvis. It organizes a wider and structured vision of space, enhances the perception of distances and environment that surrounds the user. It develops physical self-sufficiency and psychological independence by experimenting the possibility of acting independently without the help of other people.

Grillo compensates the user lack of balance thanks to its ergonomic support and wide base, it allows the user to concentrate his/her residual abilities on walking.

The wide range of components and supports offers maximum versatility. All accessories are interchangeable quickly and without the use of tools to suit different clinical pictures.



Product registered FDA (USA)

Product marked 



- **LIGHT ALUMINUM FRAME**, foldable for an easy transportation with **PATENTED HEIGHT ADJUSTMENT SYSTEM**.
- **PELVIC SUPPORT and THORACIC SUPPORT**, padded, adjustable in height, width and circumference. They allow an easy access. They hold the user just like two HANDS giving him/her a sense of security and stability allowing the user to concentrate on walking without any fear. The Thoracic support is adjustable in tilt and backward-forward.
- **ERGONOMIC HARNESS**, comfortable, adjustable and easily removable for washing. It supports the user in case of necessity without generating painful pressures.
- **HANDLEBAR**, removable and adjustable in height, it is used as a reference and as an armrest.

- **ARM SUPPORTS**, adjustable in height, tilting, rotation and depth.
- Exclusive **ORMESA DESIGN WHEELS**, light and puncture proof : 2 front swivelling, Ø15 cm (sizes mini and Small) and Ø17.5 cm (sizes Medium and Large) and 2 rear fixed, Ø20 cm (sizes mini and Small) and Ø25 cm (sizes Medium and Large).
- **DIRECTIONAL LOCK** on front wheels, **ADJUSTABLE DRAG, ANTI-ROLLBACK DEVICE** and **WHEELS STOP** on rear wheels, easy to operate by the red pedal
- **SEMIELASTIC BUMPERS**.



Attention!
Read carefully the user and maintenance manual

FRONT DRIVE

Inner circumference
pelvic and thoracic support
from 60 to 86 cm
from 23,6 to 33,8 in



WEIGHT OF THE AID 19 kg / 41,8 lb
MAX LOAD 45 kg / 99,2 lb

REAR DRIVE


Inner circumference
pelvic and thoracic support
from 60 to 86 cm
from 23,6 to 33,8 in



WEIGHT OF THE AID 18,5 kg / 40,7 lb
MAX LOAD 45 kg / 99,2 lb

USER DIMENSIONS:

USER HEIGHT (C):
from 90 to 135 cm / from 35,4 to 53,1 in

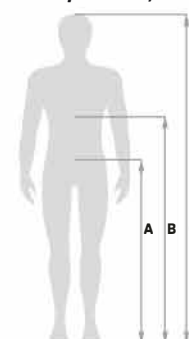

IMPORTANT:
Make sure that measures A and B indicated in the figure do not exceed the suggested range indicated for the supports.

Dealer name

Date of order

Order ref./name of the User

ORDER FORM

GRILLO



Front Drive

Rear Drive

CHOOSE THE CONFIGURATION >	PT	PA	P	PT	PA	P
ADDITIONAL COMPONENTS ▾ 925, 924: you can choose other size	 Required: > thoracic support > pelvic support > ergonomic harness > handlebar	 Required: > pelvic support > ergonomic harness > arm supports	 Required: > pelvic support > ergonomic harness > handlebar	 Required: > thoracic support > pelvic support > ergonomic harness > handles	 Required: > pelvic support > ergonomic harness > arm supports	 Required: > pelvic support > ergonomic harness > handles
865 Multiaadjustable headrest	NA	NA	NA	☐ available in combination with thoracic support 925	☐ available in combination with thoracic support 925	☐ available in combination with thoracic support 925
925 Thoracic support	●	☐	☐	●	☐	☐
924 Pelvic support	●	●	●	●	●	●
930 Ergonomic harness	●	●	●	●	●	●
926 Adjustable handlebar	●	☐	●	NA	NA	NA
809 Arm supports	☐	●	☐	☐ with joints	●	☐ with joints
943 Armstraps	☐ available in combination with arm supports 809	☐	☐ available in combination with arm supports 809	☐ available in combination with arm supports 809	☐	☐ available in combination with arm supports 809
939 Ergonomic handles	☐	●	☐	☐ with joints	●	☐ with joints
927 Removable knobs	NA	NA	NA	●	☐	●
919 Lumbar thrust support	NA	NA	NA	☐	☐	☐
918 Folding seat	NA	NA	NA	☐	☐	☐
815 Proximal abductor	☐	☐	☐	☐	☐	☐
946 Ergonomic saddle	☐	☐	☐	☐	☐	☐
890-C Leg divider with thighloops	☐	☐	☐	☐	☐	☐
890-SC Distal abductor	☐	☐	☐	☐	☐	☐
810 Weighting bars	☐	☐	☐	☐	☐	☐
923 Additional assistant push handle	☐	☐	☐	☐	☐	☐

MKT.F110 | P25001 | USA | January 2026 | R02 ● > REQUIRED IN THE CHOSEN CONFIGURATION ☐ > ACCESSORY NA > NOT AVAILABLE